



ASAM
Health is wealth. leave it to us

BATIS

Anterior Cervical Discectomy Fusion
Cervical PEEK Spacer System




Introduction

The material used is biocompatible PEEK, PEEK stands for 10993-1, Polyether Ether Ketone. PEEK polymer complies with ISO USP Class VI and ASTM F2026 for use as a medical implant material. The use of PEEK as an orthopedic device material enjoys increased popularity in recent years due to the material's unique combination of characteristics. Its properties include radiolucency, high mechanical strength, biocompatibility and compatibility with standard sterilization methods. The intrinsic radioscopic transparency of the material provides permeability on X-rays and CT scans, allowing to visualize bone growth adjacent to the implant. This allows quick and simple assessment of the bone structure and progress towards bone fusion.


Features

Anatomical shape and serrated profile for an exact implant fit and high primary stability
Optimized ratio between contact area and opening
Option of filling with bone or bone substitute to enhance bone bridging
Adequate range of sizes providing the right implant to fit the patient
Integrated blade locking mechanism provides a more reliable hold on to the endplates
The blade is designed to support primary and secondary stability to the cage.
Strong grip and easy remove

Integrated blade locking mechanism

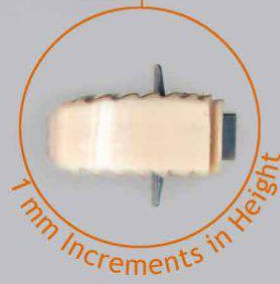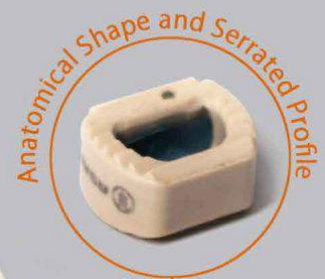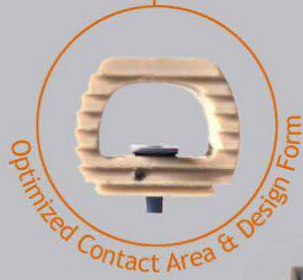


BATIS PEEK Cage-ordering information

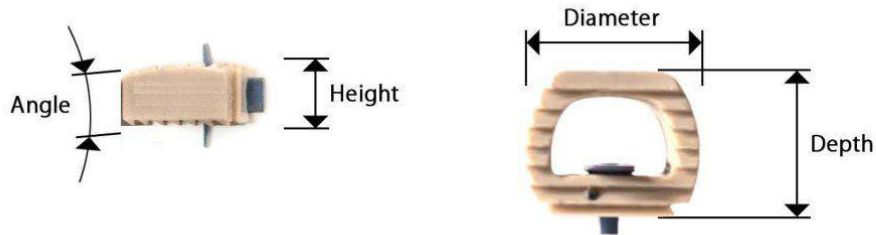

Ref No.	Height	Diameter	Depth
11-141104	4	14	11.5
11-141105	5	14	11.5
11-141106	6	14	11.5
11-141107	7	14	11.5
11-141108	8	14	11.5
11-141109	9	14	11.5
11-161304	4	16	13.5
11-161305	5	16	13.5
11-161306	6	16	13.5
11-161307	7	16	13.5
11-161308	8	16	13.5
11-161309	9	16	13.5

Please visit [www. asammedical.com](http://www.asammedical.com) for full information on indications for use, contraindications, warnings, precautions, adverse reactions and sterilization.



 No. 104, Between 5 & 7th Talash
St, Kashmar Industrial Estate, Kashmar,
Khorasan Razavi Province, IRAN

www.asammedical.com
info@asammedical.com
Document No: ASM-PC-QA-001
Date: 2023.01.07